

CARDIOTROPHIN-LIKE CYTOKINE FACTOR 1 (CLCF1) (H) ELISA KIT

FOR THE QUANTITATIVE DETERMINATION OF
HUMAN CLCF1 CONCENTRATIONS IN CELL
CULTURE SUPERNATES, SERUM, AND PLASMA



ALWAYS REFER TO LOT SPECIFIC PROTOCOL
PROVIDED WITH EACH KIT FOR
INSTRUCTIONS. PROTOCOL MUST BE
READ BEFORE USING THIS PRODUCT.

FOR RESEARCH USE ONLY. NOT FOR USE IN
DIAGNOSTIC PROCEDURES.

PRODUCT INFORMATION:

ELISA NAME	CLCF1 (H) ELISA
Catalog No.	SK000158-01
Formulation	96 T
Lot No.	
Standard range	3.125 - 200 ng/mL
Sensitivity	50 pg/mL
Sample Volume	100 µL
Dilution Factor	Optimal dilutions should be determined by each laboratory for each application
Sample Type	Serum, EDTA Plasma, Cell Culture Supernates
Specificity	Human CLCF1
Calibration	Human CLCF1 recombinant
Intra-assay Precision	4 - 6%
Inter-assay Precision	8 - 10%
Storage	2 – 8° C
This kit contains sufficient materials to run 35 samples duplicated provided that assay is run according to protocol.	

Order Contact:

AVISCEIRA BIOSCIENCE, INC.
2348 Walsh Ave., Suite C
Santa Clara, CA 95051
USA

Tel: (408) 982 0300

Fax: (408) 982 0301

Email: Sales@AvisceraBioscience.com

Info@AvisceraBioscience.com

www.AvisceraBioscience.com

DESCRIPTION

This Human CLCF1 ELISA Kit contains the necessary components required for the quantitative measurement of recombinant and/or natural human CLCF1 from cell culture supernates, serum and plasma in a sandwich ELISA format.

This immunoassay contains recombinant human CLCF1 and antibodies raised against this protein. Results from this immunoassay have shown to accurately quantify recombinant and natural CLCF1 samples.

ASSAY OVERVIEW

This assay employs the quantitative sandwich ELISA format. The plate is pre-coated with an antibody specific for human CLCF1. The capture antibody can bind to the human CLCF1 in the standard and samples. After washing the plate of any unbound substances, an antibody against human CLCF1 is added to the wells. After another washing of the plate, Anti Rabbit IgG-HRP Conjugate is added. After the last wash to remove any unbound enzyme, a substrate solution (TMB) is added to the wells and color develops in direct proportion to the amount of human CLCF1 bound in the standard solutions or samples. A standard curve can be established and sample values can be read off the standard curve.

PROCEDURAL LIMITATIONS

_FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC PROCEDURES.

_This ELISA kit should not be used beyond the expiration date on the kit label.

_Do not mix reagents with those from other lots or sources.

_It is important that the Dilution Buffer selected for the standard curve be consistent with the samples being assayed.

_Each laboratory must determine the optimal dilution factors for the samples being assayed with a pretest. If samples generate values that are not within the dynamic range of the standard curve, further concentrate or dilute the samples as required with Dilution Buffer and repeat the assay.

_Any modifications in buffers, pipetting technique, washing technique, incubation time or temperature, as well as kit age can cause a change in signal.

_Not all interfering factors have been tested in the immunoassay, therefore the possibility of interference cannot be excluded.

COMPONENTS PROVIDED

DESCRIPTION	CODE	QUANTITY
CLCF1 Microplate – 96 well microplate coated with an antibody specific for human CLCF1.	158-01-01	1 plate
CLCF1 Standard – 1000 ng/vial of lyophilized recombinant human CLCF1.	158-01-02	1 vial
Detection Antibody Concentrate – 1.05 mL/vial of 10-fold concentrate of lyophilized antibody against human CLCF1.	158-01-03	1 vial
Positive Control – one vial of lyophilized recombinant human CLCF1.	158-01-04	1 vial
Anti Rabbit IgG-HRP Conjugate – 120 µL/vial of 100-fold concentrated solution of Anti Rabbit IgG-HRP conjugate.	ARIGHRP	1 vial
Dilution Buffer – 60 mL of buffered solution with preservative.	DB06	1 bottle
Antibody and HRP Diluent Solution – 30 mL of buffered solution with preservative.	DB08	1 bottle
Wash Buffer – 50 mL of 10-fold concentrated buffered surfactant with preservative.	WB01	1 bottle
Substrate Solution – 11 mL of substrate solution.	SS01	1 bottle
Stop Solution – 11 mL of 0.5M HCl.	S-STOP	1 bottle
Plate Sealer	EAPS	1 piece
Plastic Pouch	P01	1 piece

STORAGE

Unopened Kit: Store at 2 – 8° C for up to 6 months. For longer storage, unopened Standard, Positive Control and Detection Antibody Concentrate should be stored at -20° C or -70° C. Do not use kit past expiration date.

Opened / Reconstituted Reagents: Reconstituted Standard (stock) solution and Detection Antibody concentrated solution **SHOULD BE STORED** at -20° C or -70° C for up to one month. ARIGHRP Conjugate 100-fold concentrated solution and TMB Substrate Solution can be stored at 2 – 8° C for up to 6 months (**DO NOT FREEZE** and **PROTECT FROM LIGHT**). All

other components may be stored at 2 – 8° C for up to 6 months.

Microplate Wells: Return unused strips to the plastic pouch with the desiccant pack. Microplate may be stored for up to 6 months at 2 – 8° C.

ADDITIONAL MATERIALS REQUIRED

- Microplate reader capable of absorbance measurement at 450 nm.
- Microplate shaker (250 – 300 rpm).
- Microplate washer or manifold dispenser.
- 100 mL and 500 mL graduated cylinders.
- Multi-channel Pipette, Pipettes and pipette tips.
- Deionized or distilled water.

PRECAUTION

This kit should be handled by those persons who have been trained in and can follow the principles of good laboratory practice. Wear protective clothing such as laboratory overalls, safety glasses and gloves. Care should be taken while handling solutions in this kit to avoid contact with skin or eyes, especially with the stop solution because it contains diluted hydrochloric acid. Wash immediately with water in case of contact on skin or eyes.

SAMPLE COLLECTION AND STORAGE

Cell Culture Supernates – Centrifuge and assay immediately or aliquot and store samples at ≤ -20° C. Avoid repeated freeze-thaw cycles.

Serum – Use a serum separator tube (SST). Allow blood to clot for 30 minutes. Centrifuge at 1000 x g for 15 minutes and collect serum. Assay samples immediately or aliquot and store at ≤ -20° C. Avoid repeated freeze-thaw cycles.

Plasma – Collect plasma using EDTA, heparin, or citrate as an anticoagulant. Centrifuge at 1000 x g for 15 minutes and collect plasma. Assay samples immediately or aliquot and store at ≤ -20° C. Avoid repeated freeze-thaw cycles.

Optional: Use Aprotinin (enzyme inhibitor) for ALL sample collection to prevent sample degradation. 0.5 TIU per mL of sample solution.

SAMPLE PREPARATION

Serum and plasma samples may not need to be diluted. **Optimal dilutions should be determined by each laboratory for each application with a pretest. Use polypropylene test tubes.**


REAGENT PREPARATION

Bring all reagents to room temperature before use.

Wash Buffer – If crystals have formed in the concentrate, warm bottle in a water bath until the crystals have completely dissolved. Dilute 50 mL of Wash Buffer Concentrate into 450 mL distilled or deionized water to make 500 mL of 1x Wash Buffer.

CLCF1 Standard – Reconstitute the CLCF1 standard with 1.0 mL of Dilution Buffer. The concentration of the reconstituted stock solution is 1000 ng/mL. Allow the stock standard to sit for at least 15 minutes with gentle agitation until completely dissolved prior to making standard dilutions (see below). Mix each tube thoroughly before the next transfer. Create a standard curve using a 2-fold serial dilution in Dilution Buffer with a high standard of **200 ng/mL** is recommended.

TUBE	STANDARD	DILUTION BUFFER	CONCENTRATION
Stock	powder	1.0 mL	1000 ng/mL
# 1	100µL of stock	400µL	200 ng/mL
# 2	250µL of 1	250µL	100 ng/mL
# 3	250µL of 2	250µL	50 ng/mL
# 4	250µL of 3	250µL	25 ng/mL
# 5	250µL of 4	250µL	12.5 ng/mL
# 6	250µL of 5	250µL	6.25 ng/mL
# 7	250µL of 6	250µL	3.125 ng/mL

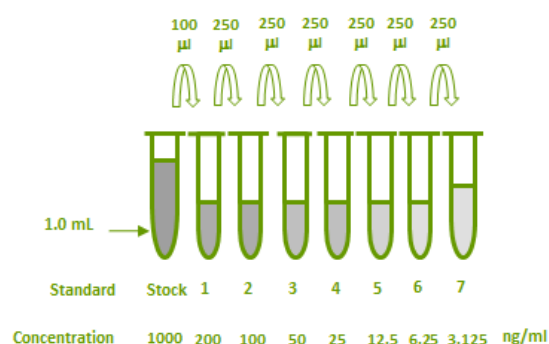

Positive Control - Reconstitute the Positive Control with 2.0 mL Dilution Buffer. **Note:** Positive Control could be reused within a few days if stored at -20° C or -70° C.

Detection Antibody - Reconstitute the Detection Antibody Concentrate with 1.05 mL of **Antibody and HRP Diluent Solution (DB08)** to produce a 10-fold

concentrated stock solution. Allow the concentrated solution to sit for at least 5 minutes until completely dissolved. Pipette 9.45 mL of Antibody and HRP Diluent Solution into a 15 mL centrifuge tube and transfer 1.05 mL of 10-fold concentrated stock solution to prepare working solution.

ARIGHRP Conjugate - Pipette 11.88 mL of **Antibody and HRP Diluent Solution (DB08)** into a 15 mL centrifuge tube and transfer 120 µL of 100-fold concentrated stock solution to prepare working solution. **Note:** 1x working solution of Anti Rabbit IgG-HRP conjugate should be used within a few days (**protect from light**). **DO NOT FREEZE.**

ELISA PROTOCOL

Bring all reagents and samples to room temperature before the start of the assay. Blank, standard dilutions, positive control and samples should be assayed in duplicate. ELISA Protocol may need further optimization.

1. Prepare all reagents, standard dilutions, positive control and samples as directed previously.
2. Remove unneeded microplate strips from the plate frame and return them to the plastic pouch with the desiccant pack.
3. Add 100 µL per well of **Dilution Buffer** to Blank wells.
4. Add 100 µL per well of **Standard dilutions, samples, or positive control**. Cover with plate sealer and incubate for 2 hours on microplate shaker at room temperature.
5. Aspirate and wash each well with 300 µL of **1x Wash Buffer** four times. After the last wash, aspirate any remaining 1x Wash Buffer, invert the plate and blot against clean paper towel(s).
6. Add 100 µL per well of **Detection Antibody working solution**. Cover with plate sealer and incubate for 2 hours on microplate shaker at room temperature.
7. Repeat the aspiration and wash as in step 5.
8. Add 100 µL per well of **ARIGHRP Conjugate working solution**. Cover with plate sealer and incubate for 45 minutes on microplate shaker at room temperature. **Protect from light.**
9. Repeat the aspiration and wash as in step 5.
10. Add 100 µL per well of **Substrate Solution**. Incubate for 10-20 minutes on microplate shaker at room temperature. **Protect from light.**

11. Add 100 µL per well of **Stop Solution**. The color in the wells should change from blue to yellow. If the color change does not appear uniform, gently tap the plate to ensure thorough mixing.
12. Read plate using a microplate reader set to 450 nm within 15 minutes.

CALCULATION OF RESULTS

Create a standard curve by plotting the log of the known concentrations of the standard dilutions (x-axis) versus the log of its corresponding O.D. (y-axis) and draw the best fit line through the points. It is recommended to use computer software capable of generating a log-log curve fit to more accurately quantify the standard dilutions.

If samples have been diluted, the concentration read from the standard curve must be multiplied by the dilution factor.

SPECIFICITY

Protein	Cross-reactivity (%)
Human CLCF1	100
Human IL-6	0
Human CNTF	0
Human CTGF	0
Human BDNF	0
Human ADRP	0









TYPICAL STANDARD CURVE

This standard curve is for demonstration only. A new standard curve should be made for each set of samples assayed.

STANDARD (NG/ML)	AVERAGE OD450 (CORRECTED)
Blank	0 (0.086)
1.563 (optional)	0.017
3.125	0.022
6.25	0.045
12.5	0.091
25	0.215
50	0.469
100	0.763
200	1.681

- Lot No.:
- Positive Control:

SUMMARY OF ASSAY PROCEDURE

PREPARE REAGENTS, SAMPLES AND STANDARD DILUTIONS

Add 100 µL of standard dilutions, samples or positive control. Cover with plate sealer and incubate 2 hours on microplate shaker at RT.

Aspirate and wash 4 times.

Add 100 µL per well of Detection Antibody working solution. Cover with plate sealer and incubate 2 hours on microplate shaker at RT.

Aspirate and wash 4 times.

Add 100 µL per well of ARIGHRP Conjugate working solution. Cover with plate sealer and incubate 45 minutes on microplate shaker at RT. Protect from light.

Aspirate and wash 4 times.

Add 100 µL per well of Substrate Solution. Incubate 10-20 min on microplate shaker at RT. Protect from light.

Add 100 µL per well of Stop Solution. Read at 450 nm within 15 minutes.