\_\_\_\_\_

# **HUMAN CXCL7/NAP-2 ELISA KIT**

FOR THE QUANTITATIVE DETERMINATION OF HUMAN CXCL7/NAP-2 CONCENTRATIONS IN SERUM AND PLASMA



ALWAYS REFER TO LOT SPECIFIC PROTCOL PROVIDED WITH EACH KIT FOR INSTRUCTIONS. PROTOCOL MUST BE READ AND CHECL ALL ITEM OF EACH KIT BEFORE USING THIS PRODUCT.

FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC PROCEDURES.

# PURCHASE INFORMATION: THIS KIT IS FOR ONE TIME USE ONLY.

ELISA NAME	HUMAN CXCL7/NAP-2 ELISA KIT
Catalog No.	SK00434-01
Lot No.	
Formulation	96 T
Standard range	125-2000 pg/mL
Sensitivity	10 pg/mL
Sample require	100 μL
Dilution Factor	Optimal dilutions should be determined by each laboratory for each application
Sample Type	Serum, Plasma
Specificity	Human CXCL7 only
Calibration	Human CXCL-7 Recombinant E. Coli derived
Intra-assay Precision	4-6%
Inter-assay Precision	4-8%
Storage	2-8°C for 1 month. See page 3 for more information
approximately 4	sufficient materials to run 0 samples duplicated provided according to protocol.

**Order Contact:** 

AVISCERA BIOSCIENCE INC. 2348 Walsh Ave., Suite C Santa Clara, CA 95051

Tel: (408) 982 0300 Fax: (408) 982 0301

Email: Sales@AvisceraBioscience.com

Info@AvisceraBioscience.com

www.AvisceraBioscience.com www.AvisceraBioscience.net

#### **INTRODUCTION**

Human CXCL7 /NAP-2 ELISA Kit is a 3.5 - 4.5 hour solid phase ELISA designed to measure Human CXCL7 in serum and plasma. It contains recombinant human CXCL7 and antibodies raised against this protein. It has been shown to accurately quantify recombinant human CXCL7. Results obtained with naturally occurring CXCL7 samples showed linear curves that were parallel to the standard curves obtained using the kit standards. These results indicate that the immunoassay kit can be used to determine relative mass values for natural human CXCL7.

#### PRINCIPLE OF THE ASSAY

This assay employs the quantitative sandwich enzyme immunoassay technique. A monoclonal antibody specific for CXCL7 has been pre-coated onto a microplate. Standards and samples are pipetted into the wells and any CXCL7 present is bound by the immobilized antibody. After washing away any unbound substances, a biotinylated polyclonal antibody specific for CXCL7 is added to the wells. Following a wash to remove any unbound antibody reagent, HRP link Streptavidin is added to the wells. After washing away any unbound enzyme, a substrate solution is added to the wells and color develops in proportion to the amount of CXCL7 bound in the initial step. The color development is stopped and the intensity of the color is measured.

# LIMITATIONS OF THE PROCEDURE

- \_ FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC PROCEDURES.
- \_ The kit should not be used beyond the expiration date on the kit label.
- \_ Do not mix or substitute reagents with those from other lots or sources.
- \_ It is important that the Dilution Buffer selected for the standard curve be consistent with the samples being assayed.
- \_ If samples generate values higher than the highest standard, dilute the samples with Dilution Buffer and repeat the assay.
- \_ Any variation in standard diluent, operator, pipetting technique, washing technique, incubation time or temperature, and kit age can cause variation in binding.
- \_ This assay is designed to eliminate interference by soluble receptors, binding proteins, and other factors present in biological samples. Until all

factors have been tested in the immunoassay, the possibility of interference cannot be excluded.

# **MATERIALS PROVIDED**

DESCRIPTION	CODE	QUANTITY
CXCL7 Microplate - 96 well	434-	1 plate
polystyrene microplate (12 strips of 8 wells) coated with	01 01	·
a purified monoclonal	01-01	
antibody against CXCL7.		
CXCL7 Standard – refer to	434-	1 vial
lot of recombinant human		1 viui
CXCL7 in a buffered protein base with preservatives;	01-02	
lyophilized.		
Detection Antibody – 105	42.4	ا د ا د
μL/vial, 100-fold	434-	1 vial
Concentrate of Biotinylated	01-03	
polyclonal antibody against		
CXCL7 with preservatives; lyophilized.		
Positive Control – one vial		
of recombinant human	434-	1 vial
CXCL7, lyophilized	01-04	
Streptavidin-HRP Conjugate	SAHRP	1 vial
– 120 μL/vial, 100-fold	SATINE	1 Viai
concentrated solution of		
Streptavidin conjugate to HRP		
Dilution Buffer - 45 mL		
/vial of buffered protein	DB02	2 vials
based solution with		
preservatives		
HRP Diluent Solution -12 mL	DB08C	1 vial
/vial of buffered protein based solution with		
preservatives		
Wash Buffer – 50 mL/vial,		
10-fold concentrated	WB01	1 vial
buffered surfactant, with		
preservative.		
TMB Substrate Solution –	TMB03	1 vial
11 mL/vial of TMB substrate solution		
Stop Solution - 11 mL/vial		
of 0.5M HCL	S-STOP	1 vial
Plate Sealer	EAPS	1 piece
	LAFO	T biece

#### **STORAGE**

**Unopened Kit:** Store at 2 - 8° C for up to 1 month. For longer storage for up to 10 months, unopened Standard, Positive Control, Detection Antibody Concentrate, Dilution Buffer and HRP Diluent Solution should be stored at -20 °C. Do not use kit past expiration date.

Streptavidin- HRP Conjugate 100-fold Concentrate and other components may be stored at 2 - 8°C for up to 10 months.

Reconstituted Standard stock solution (8ng/mL) and reconstituted Detection Antibody Stock (100-fold concentrated) may stored at -20  $^{\circ}$ C  $^{\sim}$  -70  $^{\circ}$ C for up to a few weeks.

**Microplate Wells:** Return unused wells to the plastic bag containing the desiccant pack, reseal along entire edge of zip-seal. May be stored for up to 6 months at 2 - 8°C.

# **OTHER SUPPLIES REQUIRED**

- Microplate reader capable of measuring absorbance at 450 nm, with the correction wavelength set at 540 nm or 570 nm.
- Microplate shaker (350-400rpm).
- Pipettes and pipette tips.
- Deionized or distilled water.
- Squirt bottle, manifold dispenser, or automated microplate washer.
- 100 mL and 500 mL graduated cylinders.

# PRECAUTIONS FOR USE

All reagents should be considered as potentially hazardous. The stop solution contains diluted Hydrochloric acid. Appropriate care, therefore, should be taken while handling this solution. We therefore recommend that this product be handled only by those persons who have been trained in laboratory techniques and that it is used in accordance with the principles of good laboratory practice. Wear suitable protective clothing such as laboratory overalls, safety glasses and gloves. Care should be taken to avoid contact with skin or eyes. In the case of contact with skin or eyes wash immediately with water.

# SAMPLE COLLECTION AND STORAGE

**Serum** - Use a serum separator tube (SST) and allow samples to clot for 30 minutes before centrifugation for 15 minutes at  $1000 \times g$ . Remove serum and assay immediately or aliquot and store samples at  $\leq -20$ °C. Avoid repeated freeze-thaw cycles.

Plasma - Collect plasma using EDTA, heparin, or citrate as an anticoagulant. Centrifuge for 15 minutes at 1000 x g within 30 minutes of collection. Assay immediately or aliquot and store samples at ≤-20°C. Avoid repeated freeze-thaw cycles.

#### SAMPLE PREPARATION

Serum sample may require 4000-8000 dilution. The suggested 50-fold dilution is 5  $\mu$ L plus 245  $\mu$ L of Dilution Buffer (DB02). The 2000-fold dilution is 10 $\mu$ L of 50-fold diluted sample solution plus 390  $\mu$ L of Dilution Buffer (DB02).

The final 4000-fold dilution is adding 50 μL of Dilution Buffer to all sample assay wells by a Multichannel Pipettor plus 50 μL of 2000-fold diluted sample solution to each sample well.

The final 8000-fold dilution is adding 75 μL of

Dilution Buffer to all sample assay wells by a
Multichannel Pipettor plus 25 μL of 2000-fold diluted sample solution to each sample well.


Optimal dilutions should be determined by each laboratory for each application.
Use polypropylene test tubes.

# **REAGENT PREPARATION**

Bring all reagents to room temperature before use. Wash Buffer - If crystals have formed in the concentrate, warm to room temperature and mix gently until the crystals have completely dissolved. Dilute 50 mL of Wash Buffer Concentrate into deionized or distilled water (450 mL) to prepare 500 mL of Wash Buffer.

**CXCL7 Standard** – Reconstitute standard with refer to lot of Dilution Buffer. This reconstitution produces a stock solution of 8000 pg/mL. Allow the standard to sit for a minimum of 15 minutes with gentle agitation prior to making dilutions. Pipette 300  $\mu$ L of Dilution Buffer into tubes #1 to #4. Use the stock solution to produce a 4-fold dilution series (below). Mix each tube thoroughly before the next transfer. The 2000 pg/mL standard serves as the high standard. The Dilution Buffer serves as the zero standard (0 pg/mL). The optional set 8000pg/mL as the highest standard

TUBE	STANDARD	DILUTION BUFFER	CONCENTRATION
stock	powder	Refer to lot	8000 pg/mL
#1	100 μL of stock	300μL	2000 pg/mL
# 2	100μL of 1	300μL	500 pg/mL
# 3	100μL of 2	300μL	125 pg/mL
# 4 (optional)	100μL of 3	300μL	31.25 pg/mL

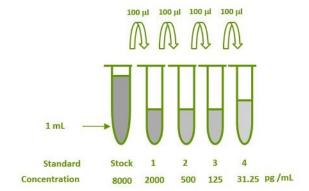

**Detection Antibody** - Reconstitute the **Detection Antibody Concentrate** with 105  $\mu$ L of Dilution Buffer to produce a 100-fold concentrated stock solution. Pipette 10.395 mL of Dilution Buffer into a 15 mL centrifuge tube and transfer 105  $\mu$ L of 100-fold concentrated stock solution to prepare working solution.

Streptavidin-HRP Conjugate - Pipette 11.88 mL of HRP Diluent Solution (DB08C) into a 15 mL centrifuge tube and transfer 120 µL of 100-fold concentrated stock solution to prepare working solution. Note: 1x working solution of Streptavidin-HRP Conjugate should be used within a few hours.

**Positive Control** - Reconstitute the **Positive Control** with refer to lot of Dilution Buffer. *Positive Control* should be prepared and used immediately. Reconstituted Positive Control CAN NOT BE REUSED.

# **ASSAY PROCEDURE**

Bring all reagents and samples to room temperature before use. It is recommended that standards be assayed in duplicate.

1. Prepare all reagents and working standards as directed in the previous sections.

- 2. Remove excess micro-plate strips from the plate frame, return them to the plastic bag containing the desiccant pack, reseal.
- 3. Add 100  $\mu$ L of **Dilution Buffer** to Blank well. **Add** 75  $\mu$ L of **Dilution Buffer** to all sample assay wells by *Multichannel Pipettor*.
- 4. Add 100 μL of Standard, or positive control per well. Add 25 μL of 2000-fold diluted sample solution to each sample assay well. (The final dilution factor for human serum samples is 8000-fold.) Cover with plate sealer. Incubate for 2 hours on micro-plate shaker at room temperature. A plate layout is provided to record standards and samples assayed.
- 5. Aspirate each well and wash, repeating the process three times for a total of four washes. Wash by filling each well with Wash Buffer (300 µL) using a squirt bottle, manifold dispenser, or autowasher. Complete removal of liquid at each step is essential to good performance. After the last wash, remove any remaining Wash Buffer by aspirating or decanting. Invert the plate and blot it against clean paper towels.
- 6. Add 100  $\mu$ L of **Detection Antibody working solution** to each well. Cover with plate sealer. Incubate for 2 hours on micro-plate shaker at room temperature.
- 7. Repeat the aspiration/wash as in step 5.
- 8. Add 100  $\mu$ L of **Streptavidin-HRP Conjugate** working solution to each well. Incubate for 45 minutes on micro-plate shaker at room temperature.
- 9. Repeat the aspiration/wash as in step 5.
- 10. Add 100  $\mu$ L of **Substrate Solution** to each well. Incubate for refer to lot at room temperature. **Protect from light.**
- 11. Add 100  $\mu$ L of **Stop Solution** to each well. The color in the wells should change from blue to yellow. If the color in the wells is green or if the color change does not appear uniform, gently tap the plate to ensure thorough mixing.
- 12. Determine the optical density of each well within 3 minutes, using a micro-plate reader set to 450 nm.

#### **CALCULATION OF RESULTS**

Average the duplicate readings for each standard, control, and sample and subtract the average zero standard optical density. Create a standard curve by reducing the data using computer software capable of generating a 4-parameter curve fit.

If samples have been diluted, the concentration read from the standard curve must be multiplied by the dilution factor.

Calculation of samples with a concentration exceeding that of standard 2000 pg/mL may result in inaccurate, low human CXCL7 levels. Such samples require further external predilution according to expected human CXCL7 values with Dilution Buffer in order to precisely quantify the actual human CXCL7 level

#### **SPECIFICITY**


Preparations of the following factors at 50 ng/mL in a mid-range rh CXCL7 control were assayed for interference. No significant cross-reactivity or interference was observed.

PROTEIN NAME	CROSS-REACTIVITY
Human CXCL7	100%
Human GRO-α	0
Human IL-8	0
Human MCP-1	0
Human SDF-1 $lpha$	0
Human SDF-1β	0

#### **TYPICAL DATA**

These standard curves are provided for demonstration only. A standard curve should be generated for each set of samples assayed.

STANDARD	AVERAGE OD450 NM		
(PG/ML)	(CORRECTED)		
BLANK	0 (refer to lot)		
125	0.069		
500	0.629		
2000	2.910		
8000 (optional)	3.643		

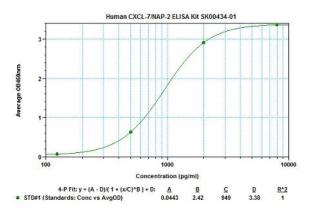

#### **SUMMARY OF ASSAY PROCEDURE**

# PREPARE REAGENTS, SAMPLES AND STANDARDS



Add 100  $\mu$ L of Dilution Buffer to blank well. Add 75  $\mu$ l of Dilution Buffer to all sample assay well by a multichannel pipettor. Add 25  $\mu$ L of 2000-fold diluted samle solution to each sample assay well. Add 100  $\mu$ L of standard, positive control to each well. Incubate 2 hours on the plate shaker at RT.



Aspirate and wash 4 times.



Add 100  $\mu$ L Detection Antibody working solution to each well. Incubate 2 hours on the plate shaker at RT.



Aspirate and wash 4 times.



Add 100 µL Streptavidin HRP conjugate working solution to each well. Incubate 45 min on the plate shaker at RT.



Aspirate and wash 4 times.



Add 100 µl Substrate Solution to each well.
Incubate 22-27 min on the bench top. **Protect from light.** 



Add 100 µl Stop Solution to each well. Read 450nm within 3 min

Use 5  $\mu$ L of human serum sample to perform 4000-fold or 8000-fold dilution of sample assay.

Dilution Factor	Samples	Dilution Buffer (DB02)
50-fold	5 μL of serum sample	245 μL
2000-fold	10 μL of 50-fold diluted sample solution	390 μL
4000-fold	50 μL of 2000- fold diluted sample solution to each sample assay well	50 μL of Dilution Buffer (DB02) to all sample assay well.
8000-fold	25 μL of 2000- fold diluted sample solution to each sample assay well	75 µL of Dilution Buffer (DB02) to all sample assay well.

The research serum samples were diluted by Dilution Buffer DB02. The linearity and recovery was assayed by SK00434-01

Sample	Dilution Factor	Assayed (pg/mL)	Final (µg/mL)	Recovery (%)
Serum F	4000 X	2034.386	8.138	100
Serum F	8000 X	990.097	7.921	97.3

The 2000-fold diluted serum sample F was assayed by SK00434-01 at 3313.697 pg/ml. That is out off the dynamic assay range 125  $^{\sim}$  2000 pg/mL. The final concentration of CXCL-7 in serum F (by 2000-fold dilution) is 6.627 µg/ml.

